

NISHIKA

Clean Away Sticky Denture Adhesive Cream Residue Quickly.

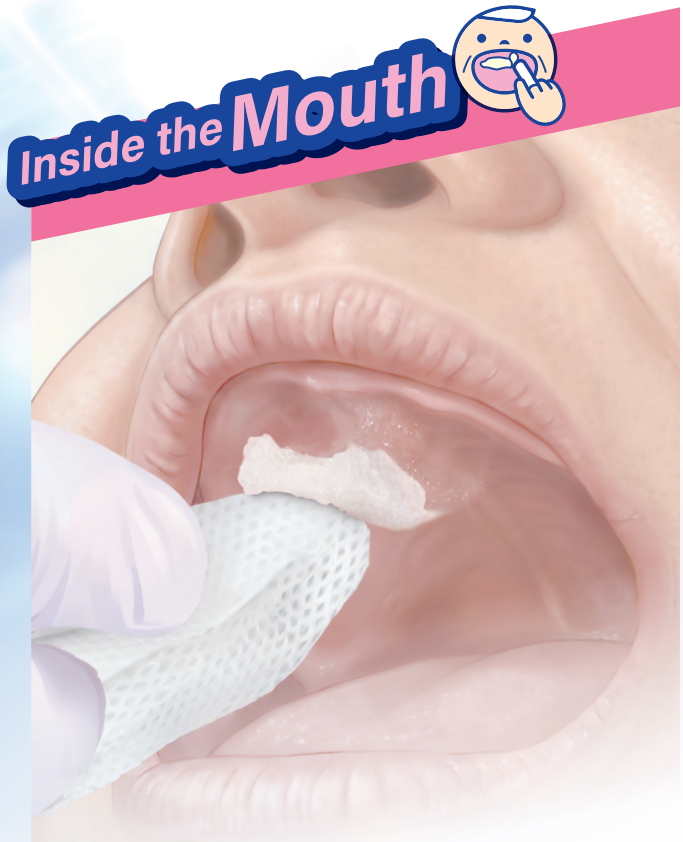

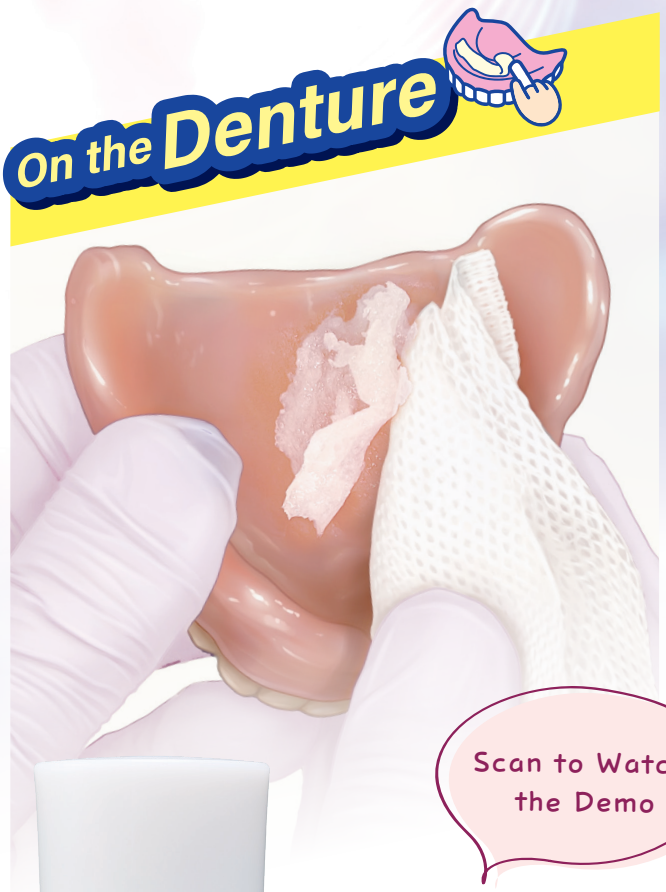

Scan to Watch the Demo

Less Effort! Less time!

Liquid quickly loses the adhesive hold, allowing denture adhesive residue to clean away easily

Denture Adhesive Remover

GRIPOFF







Bottle



Single-Dose Pouch

NISHIKA Nippon Shika Yakuhin Co.,Ltd.
4-7-24 Takezaki-cho, Shimonoseki-shi, Yamaguchi 750-0025 Japan
TEL : (81) 83-222-2221 FAX: (81) 83-222-2220
URL: https://www.nishika.co.jp/english/

A denture adhesive remover, "GRIPOFF" that easily removes cream-type denture adhesive residue clinging to the mouth and dentures, helping reduce removal time.

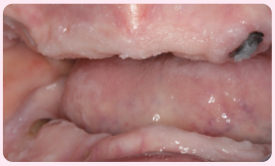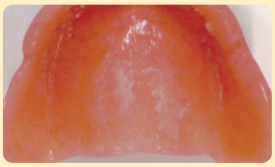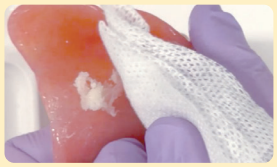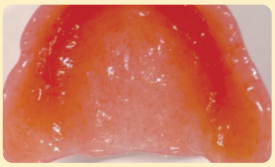	Before Removal Sticky cream-type denture adhesive residue	During Removal Lift and peel from the edge.	After Removal A refreshed appearance
Oral			
Denture			

How to use

1 Wrap the specified gauze around your finger.




2 Pour the denture adhesive remover into the tray and thoroughly soak the gauze.






3 Gently lift and peel off the denture adhesive from the edge.

Oral

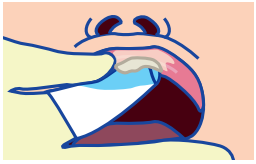


Denture

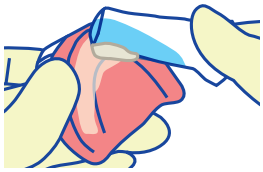

4 Repeat Steps 2–3 using a clean area of the gauze.



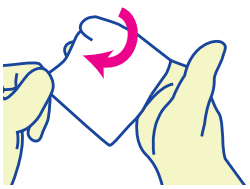

Note: If the denture adhesive gets dry and hardened, it may be difficult to remove. In this case, wet the hardened area with water first, then apply the denture adhesive remover and remove the adhesive.

- After use, rinse the mouth or the denture thoroughly with water.
- Using too little remover may result in insufficient removal.
- For small areas that cannot be cleaned with gauze, use a denturebrush or similar tool.
- A single gauze pad can be used for both the mouth and the denture.